Stay longer at home!

GRIPO ALU support pole for an easier everyday life



Functional aids for an improved quality of life See our products at hepro.no



One system – countless possibilities

GRIPO ALU a safe support in everyday life

Gripo Alu is a solid and flexible support pole that gives security when standing up and moving around in the home. Gripo Alu fits in any home, and is very easy to install. The support pole works well alone but can also be supplemented with useful accessories.

- Support when you need it
- Exploits the user's own resources
- Easy to install, move and disassemble
- Installed without screws, nails and tools
- Large selection of useful and flexible accessories
- Countless areas of application



Gripo Alu complete ART. NO. 35100



Triangel grip complete w/moss rubber and universal connector. ART. NO. 35145



• Mounted at the correct height this can also be used as headrest.

Semi hand grip

complete w/moss rubber and universal connector. ART. NO. 35159



- Simplifies rising from chair / bed
- Simplifies moving to / from wheelchail
- When using 2-3 pcs. mounted at different heights, user can train the ability to stand by using climbing effect
- By fixing the universal connection slightly loose, the user can move around the pole. Important so that user can get away from bed-chair edge

Vertical mini grip Tube diameter 22 mm. Complete w/universal connector. ART. NO. 35158







Sister assistance

complete w/2 universal connectors, hanging grip w/moss rubber and strap. ART. NO. 35148



Sister assistance

complete w/2 universal connectors, hanging grip w/moss rubber and strap. Same as picture above

Benefits of Sister assistance / Gripo Alu

- Provides support where you need it
- Provides a sense of empowerment
- Relieves nurse / helper, which can reduce strain injuries
- Activates user, which is important for among other blood
- circulation and stomach / intestinal function

• Support when moving and position-change in bed Cost saving, given that users and relatives can do more without the assistance of nursing staff

Sisterassistance tube for lateral support wo/strap and hanging grip. Complete w/2 universal connectors. ART. NO. 35196



• no loose parts

• support when standing up and moving. For example when using the toilet



Handrailtube diameter 40 mm. 75 cm length. Universal connectors ordered separately. ART. NO. 35150



- Mounted between two poles at low altitude
- Supports the user with dressing and undressing situation
- Provides user with renal function, the possibility o dressing and undressing himself

Universal connector w/joint For tube diameter 22 mm. ART. NO. 35160



Universal connector w/receiver For tube diameter 22 mm. ART. NO. 35161



Handrail tube diameter 22 mm. 75 cm length. Universal connectors ordered separately. ART. NO. 35149

Wall bracket For tube diameter 22 mm. Must be screwed into the wall. ART. NO. 35156





Combi grip Complete w/2 universal connectors. ART. NO. 35208



- Good two-hand grip for user with large hands
- Very stable. Do not twist
- Can be used for grip in two heights (climb)
- User can lean on grip while dressing and undressing


Rectangular table

41x31cm. Complete w/universal connector. Well suited as a computer-table. ART. NO. 35165



Seat Complete w/universal connector. Size: D 40 cm x W 36 cm ART. NO. 35154



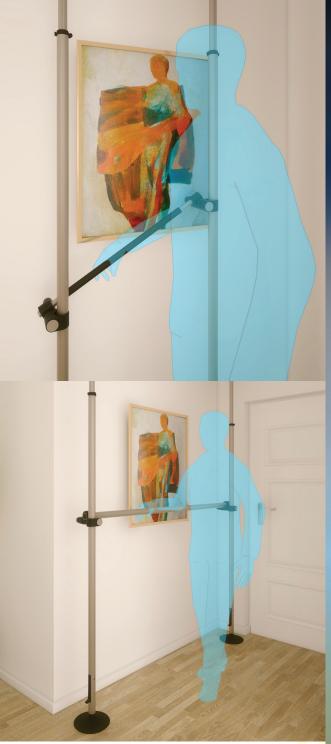



Handrail tube diameter 22 mm. 75 cm. length w/moss rubber. Universal connectors ordered separately. ART.NO 35151



- Can be installed obliquely, for example as handrails in a small stairway or where there is level difference.
- Mounted on the bathroom, by the sink or toilet between two poles.

Handrail tube diameter 40 mm.

160 cm. length. Universal connectors ordered separately. ART. NO. 35152



Mounted as walking path between two poles

- For rooms with high threshold where the user can't enter with for example wheelchair
- Walking Path for example between bed and toilet, in the corridor or such. Helps the user to manage on his own.
- Can be installed after another as long walking path between several poles.

"Because quality of life is to stay at home as long as possible"

• Proper use of Gripo Alu with different accessories can in many cases contribute so that user can stay longer at home.

Article number overview

ART.NO	ARTICLE
Gripo Alu complete	
35100	Gripo Alu, for height 210–300 cm. complete.
GRIPO ALU IS AVAILA	BLE AS SPECIAL SHORT OR LONG VERSION
35202	Gripo Alu special, for height shorter than 210 c
35105	Gripo Alu long, for height 2,6 – 3,2 mtr.

:m.



Universal connector 3!

3!

35166

5130	Universal connector	40/22 mm
5169	Universal connector	40/40 mm
5199	Universal connector	22/22 mm



Universal connector w/joint 35160 Universal connector w/joint 22 mm



Universal connector w/receiver



Wall bracket



Sister Assistance

Gripebøyle

35145

35170

35148	Sister assistance complete w/2 universal connectors,	
	hanging grip w/moss rubber and strap. (Max user weight 110 kg.)	
35133	Sister assistance w/hanging grip w/moss rubber and strap. Without	
	universal connectors. (Max user weight 110 kg.)	



Sister assistance tube for lateral support

35446	Sister assistance, tube w/clamp slider
	Without universal connectors. (Max user weight 110 kg.)
35196	Sister assistance, tube w/clamp slider and universal connectors.
	(Max user weight 110 kg.)



Sister assistance hanging grip

35132	Hanging grip w/moss rubber. Without strap.
35171	Hanging grip $w/moss$ rubber and strap.

Triangel grip w/moss rubber, complete w/universal connector. Triangel grip w/moss rubber, without universal connector.

	ART.NO	ARTICLE	
	Semi hand grip 35159 35172	Semi hand grip w/moss rubber, complete w/universal connector. Semi hand grip w/moss rubber, without universal connector.	
	Vertical mini grip		
	35158	Vertical mini grip diameeter 22 mm. Complete w/2 universal con- nectors.	
	Combi grip 35208	Combi grip complete w/2 universal connectors.	
	Handrail tube		
	35149	Handrail tube 75 cm. Diameter 22 mm.	
	35151	Handrail tube 75 cm. Diameter 22 mm. w/moss rubber	
	35150	Handrail tube 75 cm. Diameter 40 mm.	
	35152	Handrail tube 160 cm. Diameter 40 mm.	
	Cont		
	Seat 35154	Seat complete w/universal connector. (Max user weight 110 kg)	
	35173	Seat complete w/universal connector. (Max user weight no kg) Seat complete without universal connector. (Max user weight 110 kg)	
	Table	Seat complete without aniversal connector. (Max user weight no kg)	
1-	35165	Rectangular table 41x31 cm. complete w/universal connector.	
•	35155	Rectangular table 41x31 cm. complete without universal connector.	
	Top plate for sla	nting roof	
	Top plate for slanting roof		

Top plate for weak ceiling

Top plate Ellipse 70x31 cm. for weak ceiling.



Toilet roll holder35197Toilet roll holder without universal connector35198Toilet roll holder complete w/universal connector 22/22 mm.35201Toilet roll holder complete w/universal connector 40/22 mm.


Gripo Alu

Materials	Aluminium and plastic
Requirements for ceiling height	Min 210 cm - Max 300 cm
Maximum load	100 kg. (Horizontal on the middle of the pole.)
User weight	Max 150 kg
Product weight	4,7 kg. (without accessories.)

The roof construction where Gripo Alu are to be mounted must withstand a point load of 100 kg.



Notes

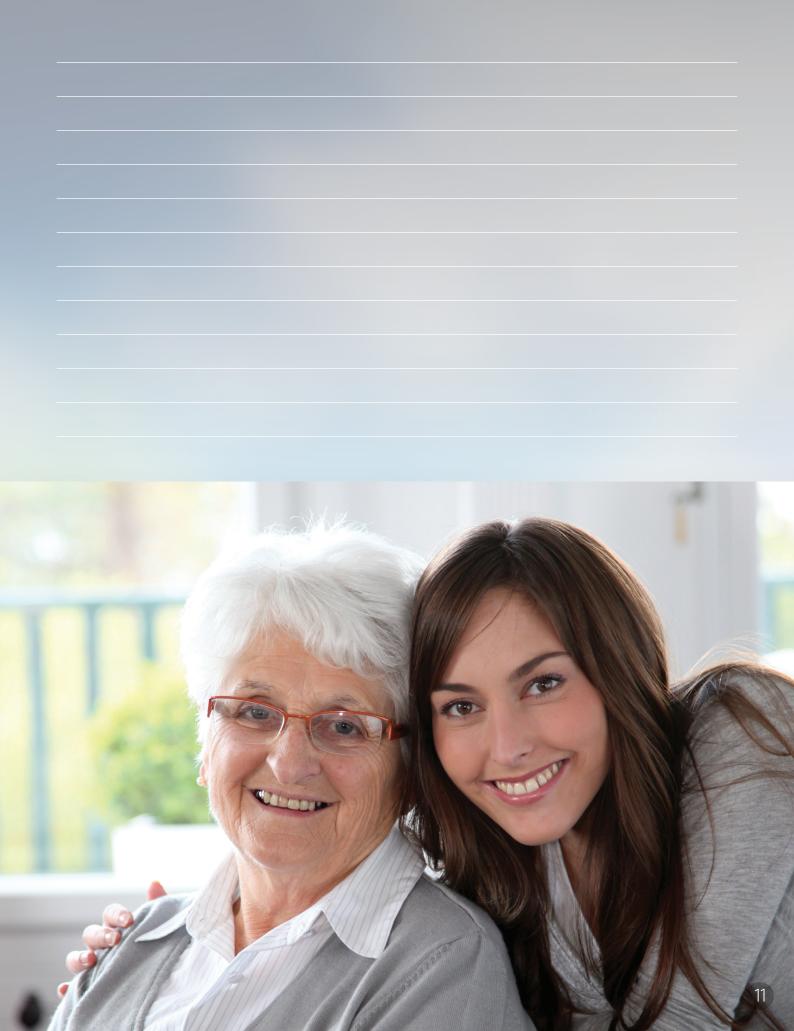

HEPRO AS Øvermoan 9, 8250 Rognan

Functional aids for an improved quality of life. See our products at hepro.no

Revisjon nr. 1 - 2015

