

AreaCare by Novis – a comprehensive range of static foam pressure relieving surfaces designed to maintain skin integrity and help prevent pressure injury, while delivering maximum comfort.

Static surfaces work by maximising the contact area for patient immersion and envelopment. The greater the area of immersion, the more the load is distributed, thereby resulting in decreased average pressure on the more vulnerable sections of the body such as heels and sacral area.

Constructed from premium materials, these products are tested to global standards to offer optimal fire retardancy and antibacterial protection. Choose from a range of specialist foam layers and features for a solution for every care environment.





As specialists in the supply of pressure care and manual transfer equipment solutions designed to help prevent pressure wound and skin integrity issues, Novis support a range of clinical healthcare products designed for hospitals, aged care facilities, rehabilitation providers and the home care community.

With a long history delivering effective pressure care management solutions, we are committed to delivering successful patient outcomes through product innovation, professional collaboration and clinical excellence.

Our product portfolio has been engineered with the goal of enhancing the comfort and safety of patients; reducing the time and concerns of care staff; and improving clinical outcomes for care facilities.

“Our goal is to provide customers with a superior choice of specialist products, while actively working to improve the quality of patient care through a program of clinical education”



Enjoy the Benefits

- Extended buying power and established infrastructure capabilities of an integrated group network.
- Specialised product knowledge and development experience, based on the latest thinking in best clinical practice for the prevention and treatment of pressure injuries.
- Comprehensive service offering – product repairs and service supported by leading warranty programs.
- Clinical expertise and industry leading training programs, including extensive understanding of clinical applications and patient outcomes.
- Fully stocked spare parts and accessories, with same-day dispatch

Learn more at [novis.com.au](https://www.novis.com.au)



AreaCare Pro

- Multi-density, three layer, zoned design for maximum comfort
- Hinged side walls which aid bed profiling and help prevent falls
- Castellated memory foam top to aid pressure redistribution and articulates with movement

CAPACITY	220 kg
SIZE	Standard



AreaCare Plus

- Multi-density, three layer, zoned design for maximum comfort
- Hinged side walls which aid bed profiling and help prevent falls
- Extra height (18cm) King Single version for greater weight capacity and support

CAPACITY	220 kg
SIZE	Standard
CAPACITY	250 kg
SIZE	King Single



AreaCare

- Multi-density, three layer design
- Memory foam top layer for optimal pressure redistribution
- Hinged side walls which aid bed profiling and help prevent falls (King Single only)


CAPACITY	200 kg
SIZE	Standard
CAPACITY	220 kg
SIZE	King Single



AreaCare Single Layer

- Economical, low risk pressure surface
- Medium density foam with inner cotton lining

CAPACITY	130 kg
SIZE	Standard

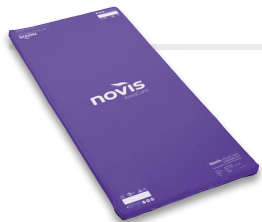

Underlay

- Provides “bottoming out” protection for alternating pressure overlays
- Minimises height of overall support system

CAPACITY	NA*
SIZE	Standard



Extension

- 200mm extension for beds with sleeping platform extension
- Memory foam top layer provides additional pressure relief under ‘at risk’ heel zone

CAPACITY	NA*
SIZE	Standard or King Single (18cm)

* Not to be used as a standalone product

Top Cover

Latex free, multi-stretch polyurethane -coated top cover for reduced shear and friction. Vapour permeable, bacteriostatically treated, fire retardant and water resistant. Double-layered waterfall skirt provides protection from moisture ingress.

Cover includes an asset management panel and turning schedule (where applicable).

Infection Control

Top cover seams are welded and/or tape sealed for maximum infection control.

AreaCare Pro, AreaCare Plus and AreaCare King Single feature a double-layered waterfall skirt to direct fluids away from mattress base.

Inner Cotton Cover

Foam covered by a removable inner cotton cover for added protection.


Hinged Base

Six hinge system to conform to profiling beds.


Base Cover

AreaCare Pro, AreaCare Plus and AreaCare King Single feature a hard-wearing, welded textured base cover material is water-resistant and features a 360° YKK zipper for easy removal.

Two carrying handles per side makes handling easy.

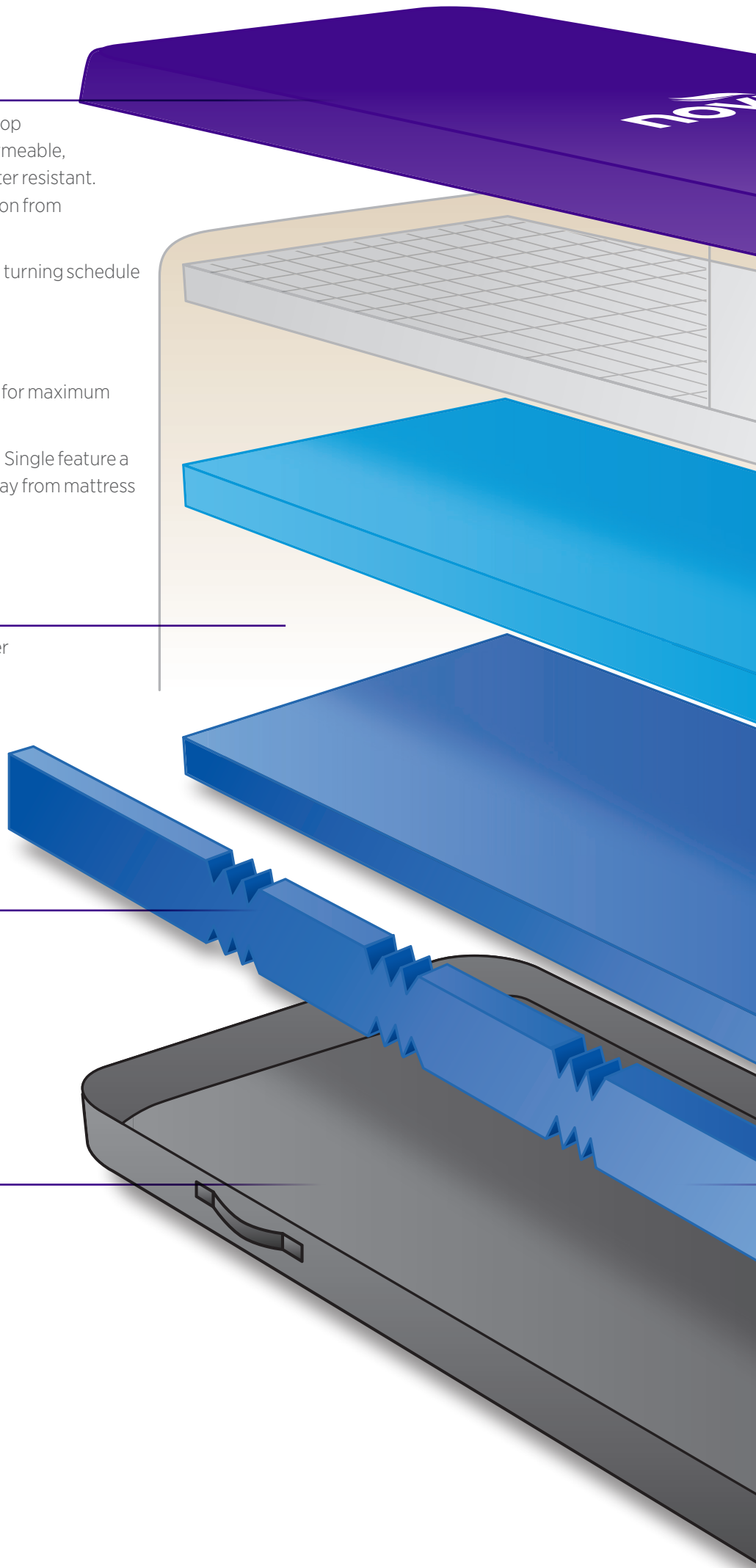

An extensive range of features to offer every patient optimum care and comfort

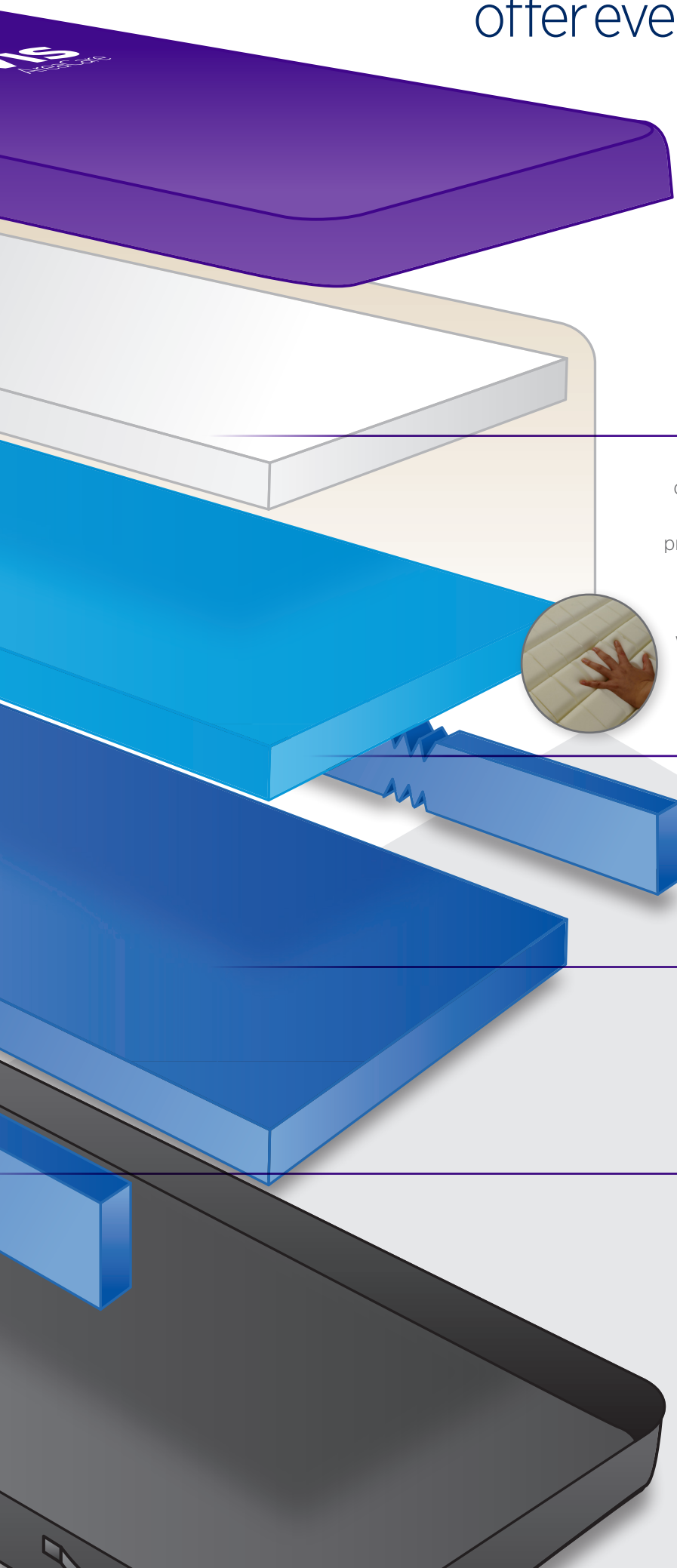

Multi-density foam layers for maximum strength and comfort – made from premium material, each foam layer is made of varying densities to maximise the redistribution of pressure.

Top Layer

High density, temperature sensitive, memory foam conforms to the body's shape; working in combination with the higher resistance support layer to disperse pressure evenly. Special immersion, slow recovery foam foot zone for the vulnerable heel area.

AreaCare Pro features a castellated memory foam top which aids pressure redistribution and articulates with movement to help reduce shear.



Support Layer

High resilience foam for maximum pressure redistribution.

Core Layer

Firm core layer and edges provide form, strength and stability for maximum safety and reduce risk of falls.

Side Bolster

Safety side bolster provide a firm surface for reduced risk of falls when sitting on edge of mattress. Includes a six hinge system to conform to profiling beds.



AreaCare Pro AreaCare Plus AreaCare Plus King AreaCare AreaCare King AreaCare Single Layer Underlay Extension Extension King

Length	1960 mm	1960 mm	1960 mm	1960 mm	1960 mm	1960 mm	1960 mm	200 mm	200 mm
Width	865 mm	865 mm	1050 mm	865 mm	1050 mm	865 mm	865 mm	865 mm	1050 mm
Height	150 mm	150 mm	180 mm	150 mm	150 mm	125 mm	50 mm	150 mm	180 mm
Layers	3	3	3	3	3	1	1	2	2
Side Walls	✓	✓	✓	✗	✓	✗	✗	✗	✗
Turn	NA	NA	NA	2 WAY	NA	4 WAY	NA	NA	NA
SWL	220 kg	220 kg	250 kg	200 kg	220 kg	130 kg	NA	NA	NA
Infection Control Seams	✓	✓	✓	✓	✓	✓	✓	✓	✓
Hinged	✓	✓	✓	✗	✓	✗	✗	✗	✗
Castellation/Convolution	✓	✗	✗	✗	✗	✗	✗	✗	✗
Waterfall Skirt	✓	✓	✓	✓	✓	✗	✗	✗	✗
Warranty	6 year body (foam) 2 year cover								
Code	BSPS3300	BSPS3000	BSPS3001	BSPS3003	BSPS3002	BSPS1100	BSPS5200	BSPS5300	BSPS5301