

GEMINO 30

walker



EXTRA SUPPORT AND STABILITY. ALSO GREAT FOR GAIT TRAINING

If you need extra support when walking, the Gemino 30 Walker is the perfect choice for you. This lightweight rollator is fitted with fully adjustable forearm supports, enabling you to walk comfortably and securely. The integrated push handles and brakes are easy to reach and very light to operate. Combining the benefits of the Gemino 30 with enhanced support and stability makes the Gemino 30 Walker also ideal for gait training.

// FOREARM SUPPORTS FOR ENHANCED STABILITY AND SAFETY

You can easily adjust the height, depth and angle of the forearm supports to suit your needs. The brakes and push handles are fully integrated – easy to reach and light to operate. For comfortable and safe walking.

// ADVANCED BRAKING OPTIONS FOR EVEN MORE SECURITY

Looking for even more security? Complement your Gemino 30 Walker with a slowdown brake or innovative SpeedControl braking system. With SpeedControl your rollator will brake automatically when you start to accelerate. You'll keep an upright posture and walk confidently and comfortably close to your rollator.



GEMINO 60 walker



OPTIMAL SUPPORT AND STABILITY – INDOORS AND OUT

Going out a lot? With the Gemini 60 Walker you get the best of both worlds: The comfortable rolling characteristics and stability of the Gemini 60 and the additional support and safety you'd expect from a 'walker'. Adjust the forearm supports and push handles to fit your needs, and enjoy more freedom indoors - and outdoors. Naturally, you can also use this lightweight rollator for extended gait training!


// ADJUSTABLE FOREARM SUPPORTS FOR ENHANCED STABILITY

The forearm supports with integrated push handles and brakes are designed for optimal support and stability, and easy use. Adjust to suit your height and enjoy a comfortable walking position.

// COMFORTABLE AND SAFE OUTDOOR USE

The sturdy frame design coupled with extra-large soft wheels allow you to go out in total comfort and safety. Cobble streets or thresholds in your home? With the Gemini 60 Walker you'll overcome any obstacle with ease.



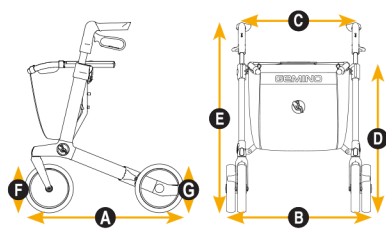

DIMENSIONS	
Maximum user weight	kg
Recommended user height	cm
Maximum load capacity in basket	kg
Total length A	cm
Total width B	cm
Total weight without basket (with basket)	kg
Width between centre of push handles C	cm
Seat height D	cm
Total height E	cm
Total height folded	cm
Total length folded	cm
Total width folded	cm
Diameter front wheels F	cm
Diameter rear wheels G	cm
PERFORMANCE	
Turning radius	cm
COLOUR OPTIONS	
Grey	
Silver grey	✓
Champagne	
Midnight blue	
Pink	
Apple green	
Black	
CERTIFICATES	
CE	✓
EN 12182	✓
ISO 11199-2	
ISO 11199-3	✓



GEMINO 30 walker

	30	30M
Maximum user weight	150	130
Recommended user height	150-200	135-170
Maximum load capacity in basket	5	5
Total length A	65	65
Total width B	61	61
Total weight without basket (with basket)	10.6 (10.9)	10.4 (10.7)
Width between centre of push handles C	25.5	25.5
Seat height D	62	55
Total height E	98 - 111.5	84 - 100
Total height folded	97	89
Total length folded	65	65
Total width folded	35.5	35.5
Diameter front wheels F	20 x 3.6	20 x 3.6
Diameter rear wheels G	20 x 3.6	20 x 3.6
PERFORMANCE		
Turning radius	84	84
COLOUR OPTIONS		
Grey		
Silver grey	✓	✓
Champagne		
Midnight blue		
Pink		
Apple green		
Black		
CERTIFICATES		
CE	✓	✓
EN 12182	✓	✓
ISO 11199-2		
ISO 11199-3	✓	✓



GEMINO 30 walker

Lightweight walker with adjustable forearm supports. The perfect companion if you need extra support and stability indoors. Also great for gait training.



Indoor use

Everyday use / gait training



150 kg
130 kg (M)



150-200 cm
135-170 cm (M)



Height-adjustable forearm supports for extra support and stability



For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.