


Pressure Relief Cushion

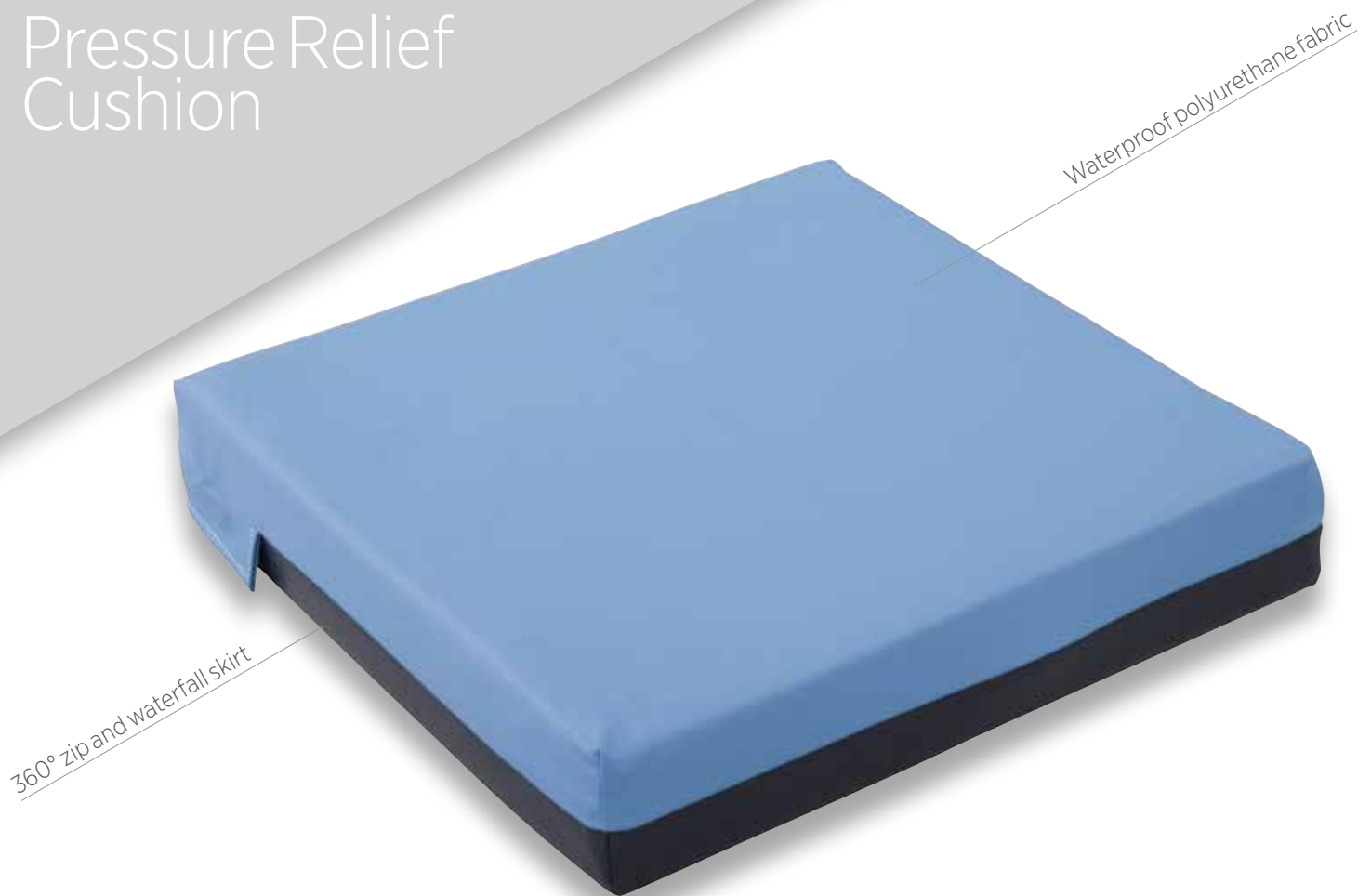

This two layer foam cushion provides improved pressure redistribution when seated.

This multi-purpose product is suitable for users with mild postural needs and comprises a high density foam support base with a super soft memory foam top layer.

Enclosed by a zippered two-way stretch waterproof polyurethane fabric cover for enhanced anti-shear and anti-friction benefits, with 360° zip and waterfall skirt to protect from moisture ingress.

Seams are stitched for strength, then welded for or maximum infection control.

o Formerly called the Comfort Cushion

CODE	BL7258
HEIGHT	75 mm
WIDTH	410 mm
DEPTH	410 mm
CAPACITY	120 kg
WARRANTY	1 year

