


EquaGel Adjustable Protector Cushion

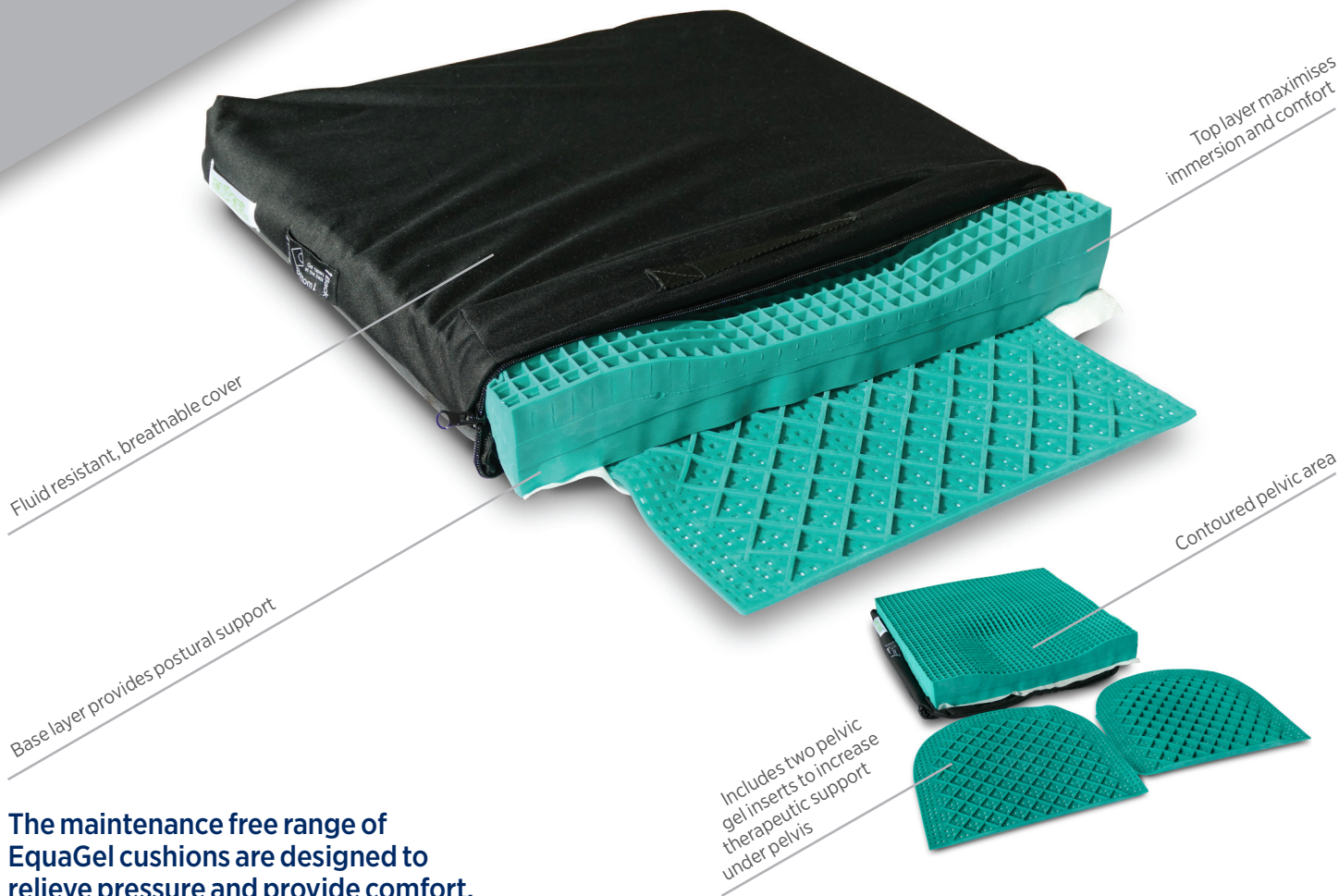

The maintenance free range of EquaGel cushions are designed to relieve pressure and provide comfort.

Features

- Includes two pelvic gel inserts to increase therapeutic support under pelvis
- Thick or thin inserts for customisable cushioning needs
- Fluid resistant, flame retardant, vapour permeable cover
- Hook and loop sections on base of cushion for attaching to a wheelchair
- 63.5mm high (without inserts) with 44.5mm of buckling gel under pressure points
- Dual layer design:
 - The top layer consists of a thinner wall gel and tight grid for superior immersion and comfort
 - The base layer consists of a thicker wall gel and wider grid which provides deep postural support
- 2 year warranty on insert
- 1 year warranty on cover
- Safe working load (SWL) 250kg

Size

Code

400 x 400 mm, 16 x 16"	8-400027-4
400 x 460 mm, 16 x 18"	8-400027-7
460 x 400 mm, 18 x 16"	8-400027-2
460 x 460 mm, 18 x 18"	8-400027-5
510 x 400 mm, 20 x 16"	8-400027-3
510 x 460 mm, 20 x 18"	8-400027-1
560 x 400 mm, 22 x 16"	8-400027-13
560 x 460 mm, 22 x 18"	8-400027-10